

CHOOSE THE RIGHT LEVEL OF PROTECTION FOR YOUR SITUATION

In today's healthcare environment, it is more important than ever for clinicians to be equipped with the proper protective apparel to help reduce the spread of Healthcare-Associated Infections (HAIs). Our Isolation Gown Portfolio provides barrier protection based on AAMI guidelines and helps reduce the risk of contamination and transmission of infectious organisms that lead to HAIs.

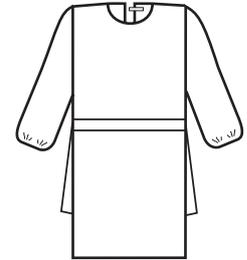
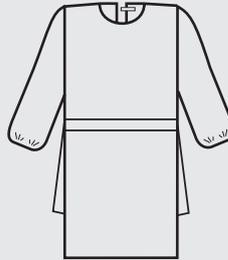


OVERVIEW PROTECTIVE GOWNS

HALYARD* ISOL

TRI-LAYER AAMI1 GOWNS

TRI-LAYER AAMI2 GOWNS



PRODUCT FEATURES

Colour	Yellow	Yellow	Yellow	Yellow	Blue	Blue
Product codes	54100	54101	69979	69988	69981	69987
Size	L	XL	L	XL	Univ	XL
Material	3 layers light-weight SMS		3 layers medium-weight SMS			
Weight	20g/m ²		29g/m ²			
ANSI/AAMI PB70 classification	1		2			
Splash resistance	Moderate		Good			
Hydrostatic pressure			≥20 cm H ₂ O			
Shelf life	3 years		3 years			

RECOMMENDED AREAS OF USE

Visitors gown						
Basic care / transport patients		✓				
Laundry/housekeeping		✓				
Patient care if C- or D- isolation					✓	
Medical/nursing acts with low risk of splashing of body fluids (ICU, dialysis, med/surg unit nursery, lab pathology, laboratories, hyperbaric...)					✓	
Medical/nursing acts with real risk of splashing of body fluids (ER, trauma, burn units, critical care units, ..)						
Surgical or lab procedures that require a non-sterile OR gown						
Procedures with a lot of fluid or liquid (instrument cleaning CSSD)						
Preparing cytostatics						

HANDLING FEATURES

Closure neck		Tape tab			Tape tab	
Closure waist		Tie			Tie	
Opening back		No			No	
Cuffs		Elastic			Elastic	

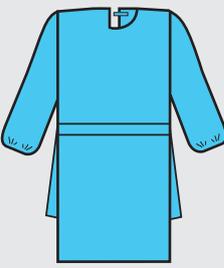
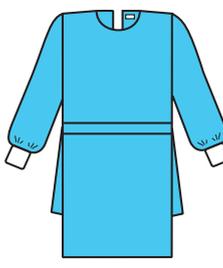
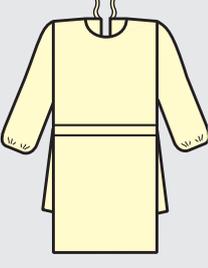
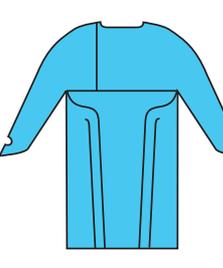
PACKAGING

100/case

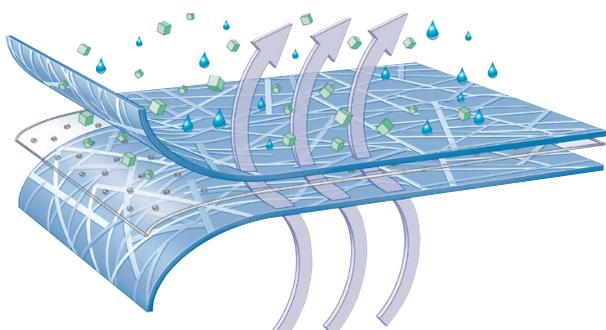
100/case

PROTECTION GOWNS

HALYARD* FILM GOWNS

TRI-LAYER AAMI3 GOWNS		PROCEDURE GOWNS		SPUNBOND GOWNS		IMPERVIOUS COMFORT GOWNS			
									
Blue	Blue	Blue		Yellow	Yellow	Blue			
54310	54311	69025	69028	69124	69125	69600	69602	69601	69603
L	XL	Univ	XL	Univ	XL	Univ	Univ	XXL	XXL
3 layers heavy-weight SMS		3 layers protective SMS		1 layer spunbond		Polymer film + spunbond			
36g/m ²		50g/m ²		20g/m ²		39g/m ²			
3									
Very good		Excellent		Weak		Impervious			
≥ 50 cm H ₂ O		84,7 cm H ₂ O				122,4 cm H ₂ O			
3 years		N/A		N/A		N/A			
				✓					
✓									
✓									
		✓						✓	✓
Tape tab		Velcro®		Tie		Pre-cut opening			
Tie		Tie		Tie		No			
No		No		No		Perforated open back			
Elastic		Knitted cuff		Elastic		Thumb hooks/Knitted cuff			
100/case		60/case		100/case		100/case			

The HALYARD* Tri-Layer Isolation Gowns are Made with Protective Tri-Layer SMS Fabric



Outer layers - Spunbond polypropylene material that is cloth-like, yet tough and abrasion-resistant

Middle layer - Intricate meltblown polypropylene microfibers that filter out many water-based liquids

AAMI Guidelines: a reliable benchmark for protection

	LEVEL 1	LEVEL 2	LEVEL 3
AAMI level / performance requirements at 4% AQL	Spray impact penetration < 4.5g	Spray impact penetration < 1.0g Hydrostatic pressure > 20cm	Spray impact penetration < 1.0g Hydrostatic pressure > 50cm
Anticipated risk of exposure to fluid	Low	Between low and moderate	Moderate
Gown characteristics (weight)	Light-weight SMS fabric	Medium-weight SMS fabric	Heavy-weight SMS fabric

To meet AAMI Guidelines, protective apparel must achieve minimum performance standards for strength, barrier performance and fluid resistance. Tested under controlled conditions, AAMI-rated gowns must demonstrate a 4% AQL level, which provides a high level of assurance that users are protected for the anticipated level of fluid contact. Each garment must be fully protective in the front, back and along the seams, and be clearly labeled to indicate the level of protection provided, from AAMI Level 1 (lowest) to AAMI Level 4 (highest).



KNOWLEDGE NETWORK* Clinical Education
Knowledgeable Customer Support
Expert Sales Force
Tools & Best Practices
Clinical Research
Commitment to Excellence

For more information, please contact your Halyard Customer Care Team:

UK & IE (toll-free)
Tel: 0800 085 82 04 - Fax: 0800 066 91 06
customerservice.uk.ie@hyh.com

South Africa (toll-free)
Tel: 0801 111 332 - Fax: 0801 111 335
customerservice.za@hyh.com

Other countries
Tel: +32 2 585 51 14 - Fax: +32 2 897 00 95
customerservice.export@hyh.com

